


Excellent Brain – Trainer Mangement System Manual

Version 1.0 – Aug. 2019

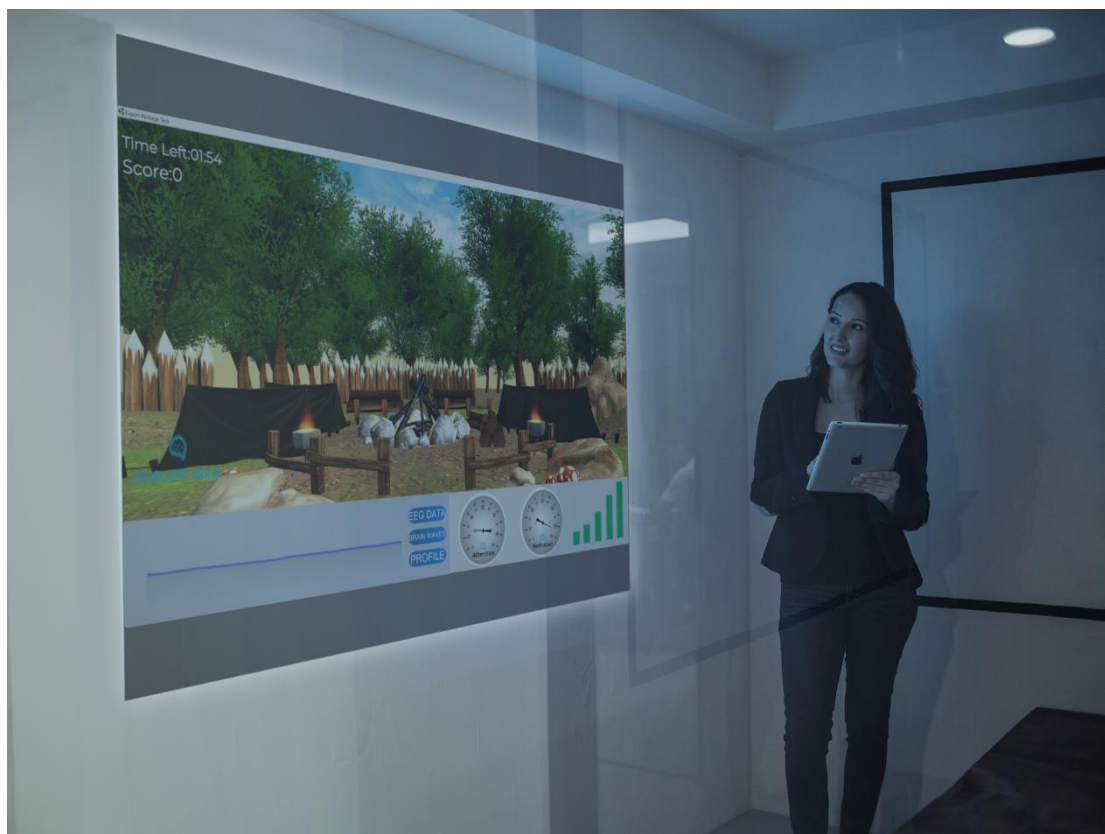


1. Reaching the website -

Go to the link –

<https://excellent-brain.co.il/Logic/View/TrainersLogic/lang/eng/Index.php>

Or from the software (version 0.95.10 and above), click *settings* in the main menu and then click the animated brain figure at the top of the window.

You will reach this website –

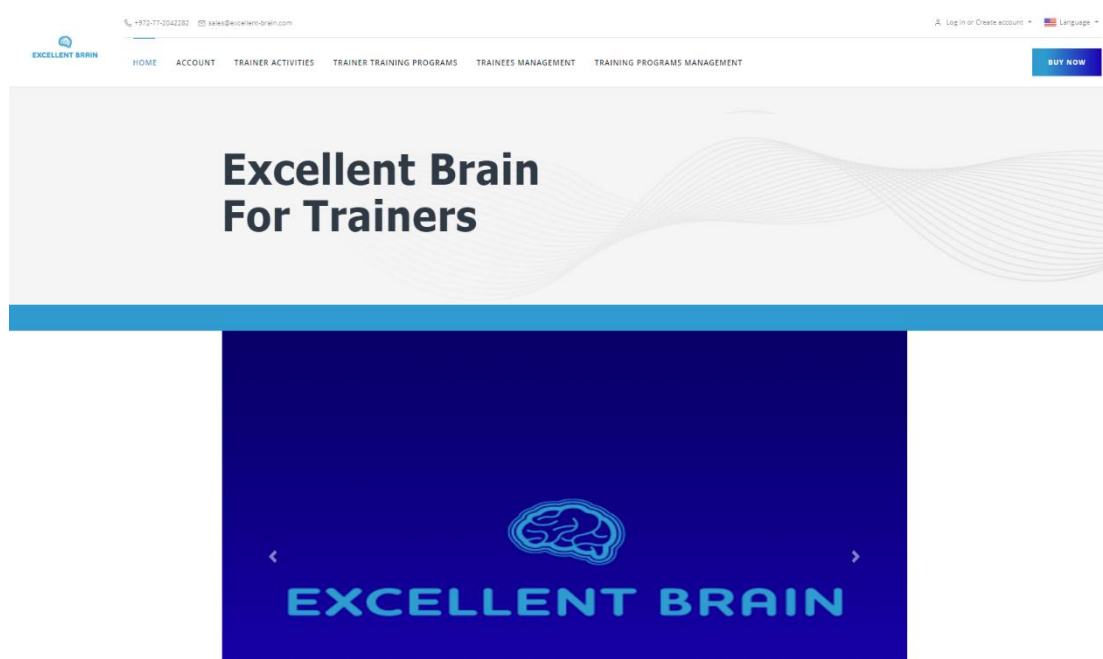


2. **Registration and Login** – Click on *Log In or Create Account* at the top of the page on the right corner, then you will reach the following page where you will register on the right side of the page.

Important! Register with your right Email address, so you will be able to get messages and follow Sessions in the future. Without the right Email address you cannot login to Excellent Brain.

Make sure to fill in all the details in the registration form and after clicking **REGISTER** wait a few seconds until you get a green message at the top of the page telling you to check your Email account for a confirmation Email.

Open your Email account and then open the Email you got from Excellent Brain and click the link you received.
That is it – Excellent Brain now recognizes you as a trainer.

Now you can go back to the website and login on the left side of the page with the Email address and password you chose before.

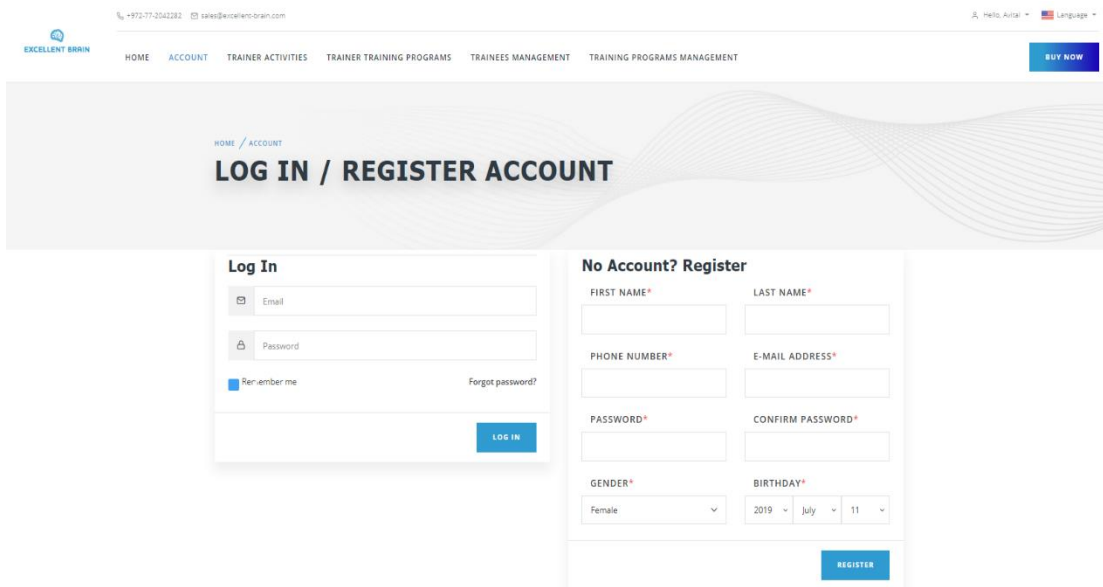

EXCELLENT BRAIN +972-77-2042282 sales@excellentbrain.com Hello, Aitor Language

HOME ACCOUNT TRAINER ACTIVITIES TRAINER TRAINING PROGRAMS TRAINEES MANAGEMENT TRAINING PROGRAMS MANAGEMENT BUY NOW

HOME / ACCOUNT

LOG IN / REGISTER ACCOUNT

Log In

Email

Password

☐ Remember me [Forgot password?](#)

LOG IN

No Account? Register

FIRST NAME* LAST NAME*

PHONE NUMBER* E-MAIL ADDRESS*


PASSWORD* CONFIRM PASSWORD*

GENDER* BIRTHDAY*

Female 2019 July 11

REGISTER

3. **Edit Personal Information** – In order to edit your personal information, click at the top of the page on *PROFILE SETTINGS*. You will reach this page -

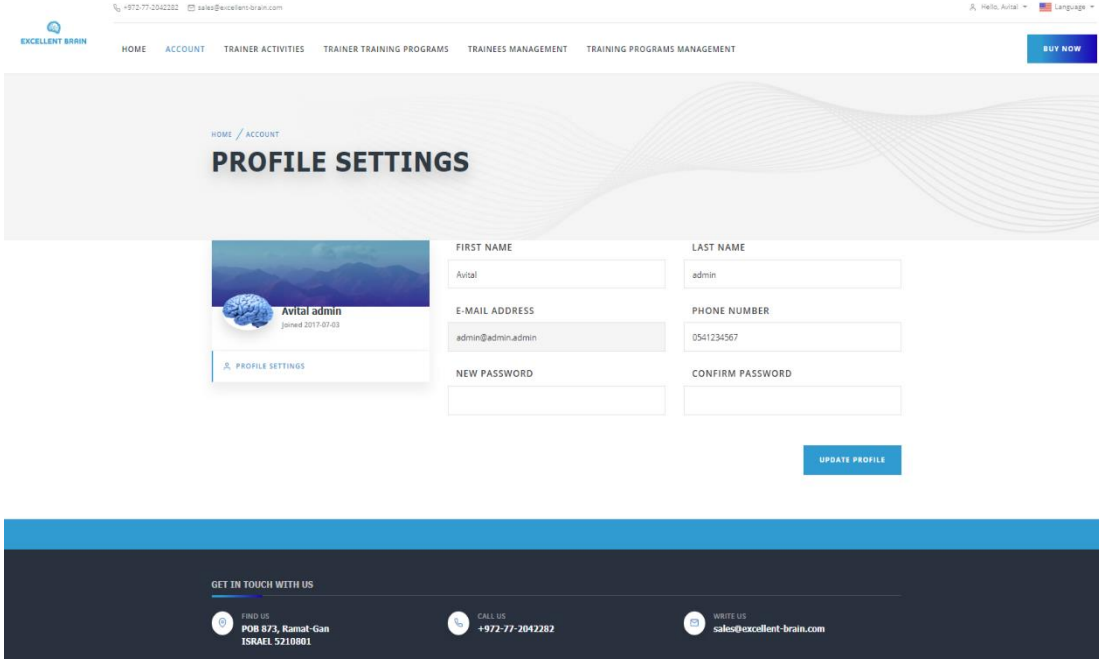

The screenshot displays the 'PROFILE SETTINGS' page. At the top, there is a navigation bar with links: HOME, ACCOUNT, TRAINER ACTIVITIES, TRAINER TRAINING PROGRAMS, TRAINEES MANAGEMENT, and TRAINING PROGRAMS MANAGEMENT. A 'BUY NOW' button is also present. Below the navigation bar, the page title 'PROFILE SETTINGS' is shown. On the left, there is a user profile card for 'Avital admin' with a profile picture and a 'PROFILE SETTINGS' link. The main content area contains a form with the following fields:


FIRST NAME Avital	LAST NAME admin
E-MAIL ADDRESS admin@admin.admin	PHONE NUMBER 0541234567
NEW PASSWORD	CONFIRM PASSWORD

An 'UPDATE PROFILE' button is located at the bottom right of the form.

At the bottom of the page, there is a footer section titled 'GET IN TOUCH WITH US' with three columns:

- FIND US:** POB 872, Ramot-Gan, ISRAEL 5210881
- CALL US:** +972-77-2042282
- WRITE US:** sales@excellent-brain.com

4. **Trainees Management** – In order to see your trainee's progress, click at the top of the page on *TRAINEES MANAGEMENT* and then you can see all of your trainees and their general data.

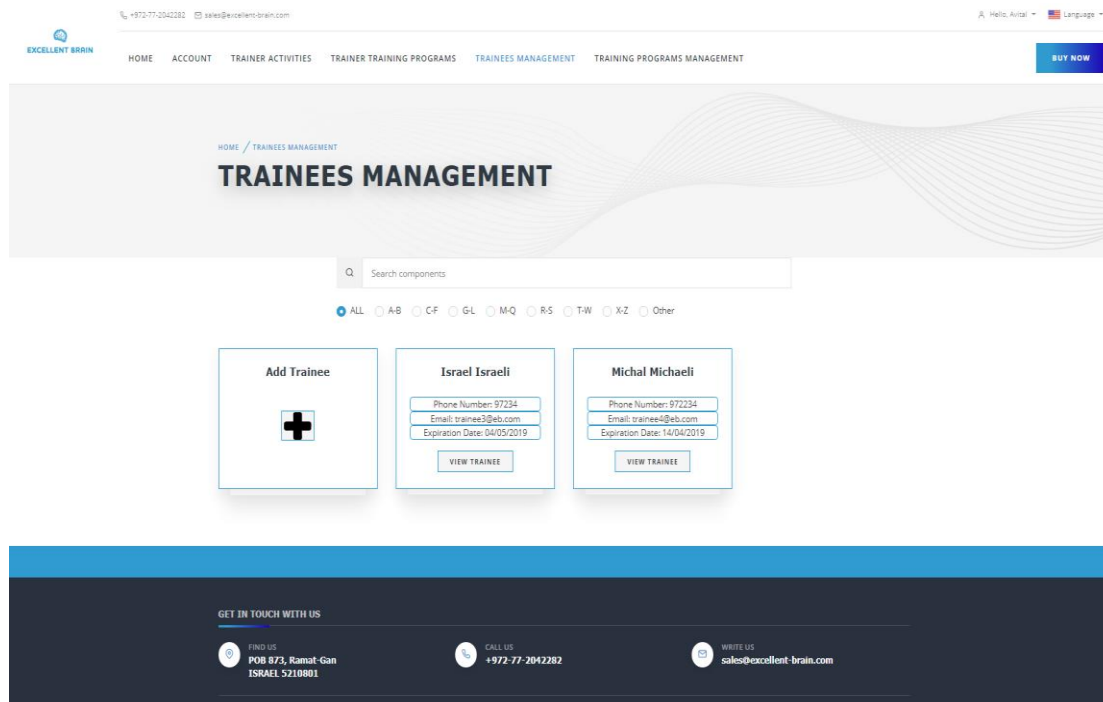

The screenshot shows the 'TRAINEES MANAGEMENT' page of the EXCELLENT BRAIN system. The page has a header with the company logo, contact information (+972-77-2042282, sales@excellent-brain.com), and a navigation menu including HOME, ACCOUNT, TRAINER ACTIVITIES, TRAINER TRAINING PROGRAMS, TRAINEES MANAGEMENT (selected), and TRAINING PROGRAMS MANAGEMENT. A 'BUY NOW' button is also present.

The main content area features a search bar for components and a filter section with radio buttons for ALL, AB, CF, GL, MQ, RS, TW, XZ, and Other. Below the filters are three cards:


- Add Trainee**: A card with a plus icon and a 'VIEW TRAINEE' button.
- Israel Israeli**: A card showing details for a trainee named Israel Israeli, including Phone Number: 97234, Email: trainee3@eb.com, and Expiration Date: 04/05/2019. It has a 'VIEW TRAINEE' button.
- Michal Michaeli**: A card showing details for a trainee named Michal Michaeli, including Phone Number: 97234, Email: trainee4@eb.com, and Expiration Date: 14/04/2019. It has a 'VIEW TRAINEE' button.

The footer section, titled 'GET IN TOUCH WITH US', provides contact information:

- FIND US**: POB 873, Ramat-Gan, ISRAEL 5210801
- CALL US**: +972-77-2042282
- WRITE US**: sales@excellent-brain.com

5. **Add New Trainees** – In order to add new trainees, click at the top of the page on *TRAINEES MANAGEMENT* and then click on the + under *Add Trainee*. Then you will move to a register page, almost identical to the one you registered from.

Important! Register with your trainee's right Email address. Without the right Email address your trainee cannot train in Excellent Brain.

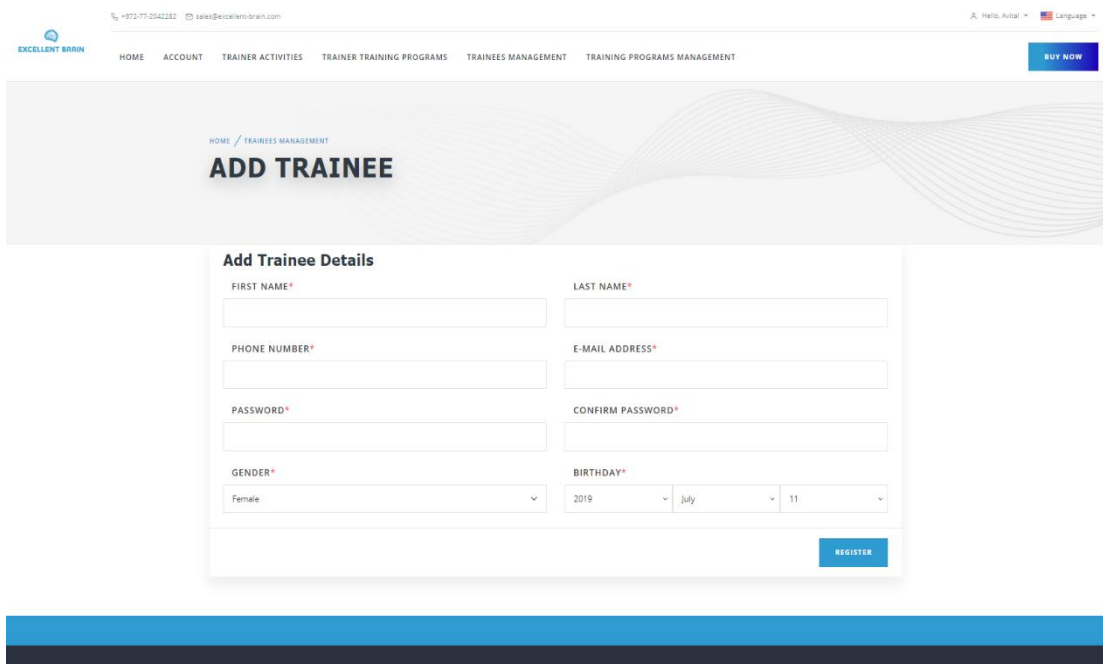


The screenshot shows the 'ADD TRAINEE' registration page. The header includes the Excellent Brain logo, contact information (+972-77-2042282, sales@excellentbrain.com), and a language selector. The navigation menu contains: HOME, ACCOUNT, TRAINER ACTIVITIES, TRAINER TRAINING PROGRAMS, TRAINEES MANAGEMENT, and TRAINING PROGRAMS MANAGEMENT. A 'BUY NOW' button is also present. The main heading is 'ADD TRAINEE'. Below it, the 'Add Trainee Details' form is displayed with the following fields:

- FIRST NAME*
- LAST NAME*
- PHONE NUMBER*
- E-MAIL ADDRESS*
- PASSWORD*
- CONFIRM PASSWORD*
- GENDER* (Dropdown menu showing 'Female')
- BIRTHDAY* (Date picker showing 2019, July, 11)

A 'REGISTER' button is located at the bottom right of the form.

6. **View All Activities Performed By the Trainee** – In order to view Activities history of a specific trainee, click at the top of the page on *TRAINEES MANAGEMENT* and then choose the trainee by clicking on *VIEW TRAINEE*. Then choose Activities option by clicking on *VIEW ACTIVITIES*.

Now you can view all Activities performed by this trainee, sorted by activity name, difficulty level and date. You can also view current score and best score of this activity in this level.

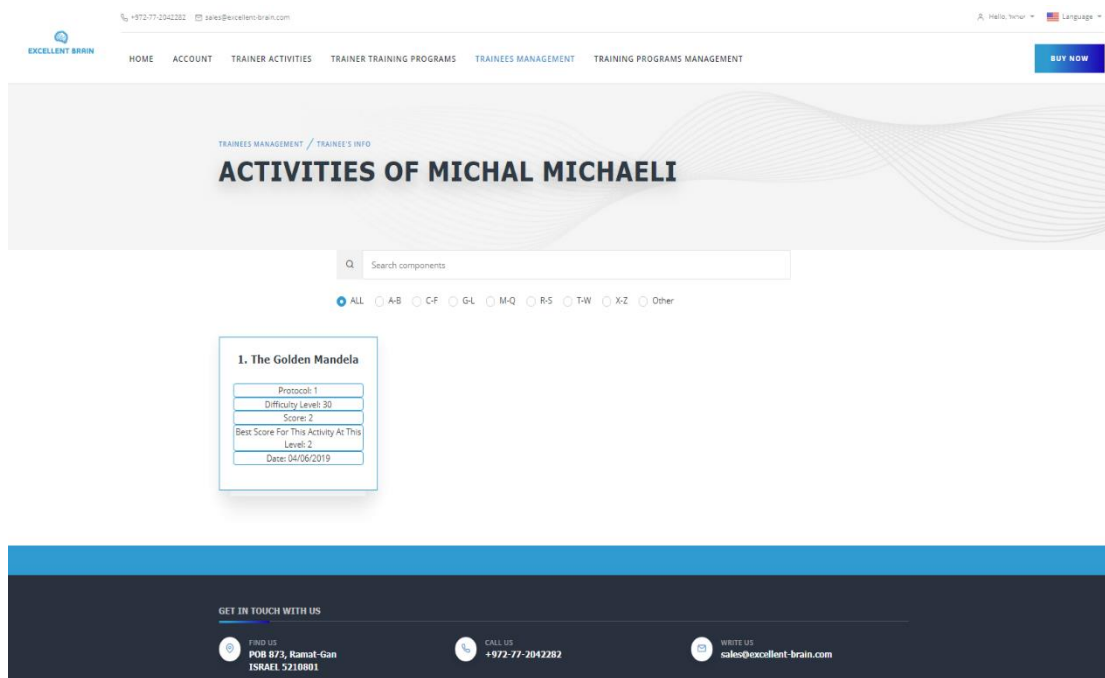

The screenshot shows the 'ACTIVITIES OF MICHAL MICHAELI' page on the Excellent Brain website. The page has a header with the company logo, contact information, and navigation links. Below the header, there is a search bar and a filter section with radio buttons for 'ALL', 'A-B', 'C-F', 'G-L', 'M-Q', 'R-S', 'T-W', 'X-Z', and 'Other'. The main content area displays a table with activity details for '1. The Golden Mandela'.


1. The Golden Mandela	
Protocol:	1
Difficulty Level:	30
Score:	2
Best Score For This Activity At This Level:	2
Date:	04/06/2019

The footer contains contact information under the heading 'GET IN TOUCH WITH US':

- FIND US:** POB 873, Ramat-Gan, ISRAEL 5210801
- CALL US:** +972-77-2042282
- WRITE US:** sales@excellent-brain.com

7. **View All Training Programs Performed By the Trainee** – In order to view Training Programs history of a specific trainee, click at the top of the page on *TRAINEES MANAGEMENT* and then choose the trainee by clicking on *VIEW TRAINEE*. Then choose activities option by clicking on *VIEW PROGRAMS*.

Now you can view all Training Programs performed by this trainee and their progress.

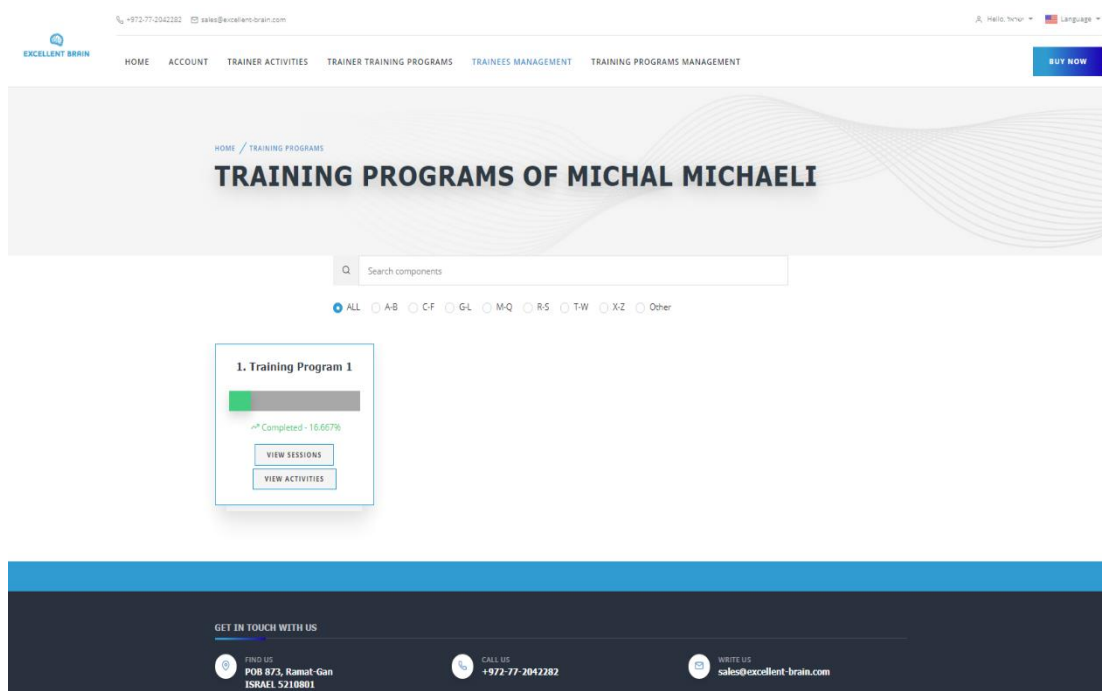


The screenshot displays the 'TRAINING PROGRAMS OF MICHAL MICHAELI' page in the Excellent Brain application. The page features a navigation bar at the top with links for HOME, ACCOUNT, TRAINER ACTIVITIES, TRAINER TRAINING PROGRAMS, TRAINEES MANAGEMENT, and TRAINING PROGRAMS MANAGEMENT. A 'BUY NOW' button is also present. Below the navigation bar, the page title 'TRAINING PROGRAMS OF MICHAL MICHAELI' is prominently displayed. A search bar labeled 'Search components' is located below the title. Underneath the search bar, there are radio buttons for filtering by category: ALL (selected), AB, CF, GL, M-Q, RS, TW, XZ, and Other. The main content area shows a list of training programs, with the first one being '1. Training Program 1'. This program has a progress bar indicating it is 'Completed - 16.667%'. Below the progress bar, there are two buttons: 'VIEW SESSIONS' and 'VIEW ACTIVITIES'.

8. **View Activities Table of Training Program Performed By the Trainee** – In order to view activities table of Training Program performed by a specific trainee, click at the top of the page on *TRAINEES MANAGEMENT* and then choose the trainee by clicking on *VIEW TRAINEE*. Then choose activities option by clicking on *VIEW PROGRAMS*.

Now you can view all Training Programs performed by this trainee and their progress.

In order to watch all Activities table, click *VIEW ACTIVITIES* and you will get this –


+972-77-204202 | sales@excellent-brain.com
Hebrew | English

[HOME](#)
[ACCOUNT](#)
[TRAINER ACTIVITIES](#)
[TRAINER TRAINING PROGRAMS](#)
[TRAINEES MANAGEMENT](#)
[TRAINING PROGRAMS MANAGEMENT](#)
[BUY NOW](#)

[TRAINEE'S INFO](#) / [TRAINEE'S TRAINING PROGRAMS](#)


ACTIVITIES IN TRAINING PROGRAM 1 - MICHAL M ICHAEI

#	Session Name	Activity Name	Place	Protocol	Difficulty Level	Is Completed
1	1	Session 1	Cupid	Attention	30%	✓
2	2	Session 1	My Earth	Attention	30%	✓
3	3	Session 1	Shark	Attention	40%	✓
4	4	Session 1	Launch the Rocket	Attention	40%	✓
5	5	Session 1	Resize the Diamond	Attention	50%	✓
6	6	Session 1	Lift the Brain	Attention	50%	✓
7	7	Session 1	Fireworks	Attention	50%	✓
8	8	Session 1	Glowing Cube	Attention	40%	✓
9	1	Session 2	Puffer fish	Attention	30%	✓
10	2	Session 2	Lift the Drone	Attention	30%	✓
11	3	Session 2	Where is the Spaceship	Attention	40%	✓
12	4	Session 2	Tunnel	Attention	40%	✓
13	5	Session 2	Cupid	Attention	50%	✓

9. **View Sessions and Sessions Reports in A Training Program Performed By the Trainee** - In order to view Sessions of a Training Program performed by a specific trainee, click at the top of the page on *TRAINEES MANAGEMENT* and then choose the trainee by clicking on *VIEW TRAINEE*. Then choose activities option by clicking on *VIEW PROGRAMS*.

Now you can view all Training Programs performed by this trainee and their progress.

To view the Sessions, click on *VIEW SESSIONS* and you will get this –

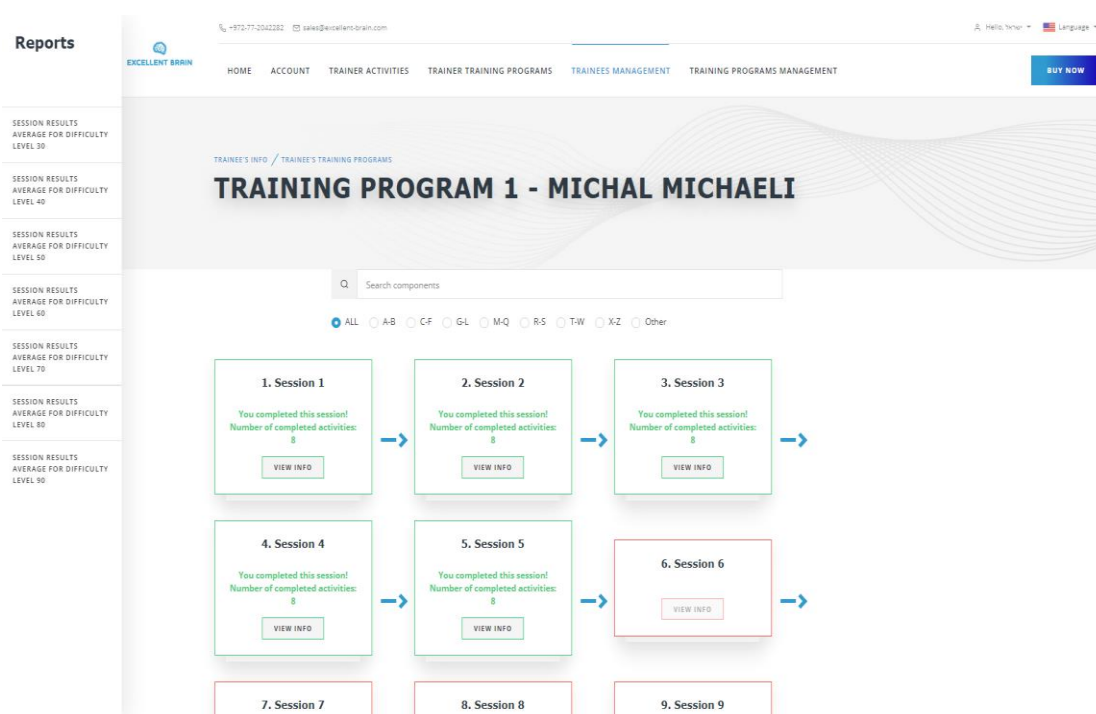


The screenshot displays the 'Reports' section of the Excellent Brain system. The main heading is 'TRAINING PROGRAM 1 - MICHAL MICHAELI'. Below this, there is a search bar and a filter menu with options: ALL, A-B, C-F, G-L, M-Q, R-S, T-W, X-Z, and Other. The sessions are arranged in a grid:

- 1. Session 1**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 2. Session 2**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 3. Session 3**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 4. Session 4**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 5. Session 5**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 6. Session 6**: You completed this session! Number of completed activities: 8. [VIEW INFO]
- 7. Session 7**: [VIEW INFO]
- 8. Session 8**: [VIEW INFO]
- 9. Session 9**: [VIEW INFO]

Arrows indicate the sequence from Session 1 to Session 9. The left sidebar shows a list of session results for various difficulty levels (30, 40, 50, 60, 70, 80, 90).

10. View a Report displaying Sessions Results Average for Difficulty Level 30 Graph (For Example) – Reach the page at section 9 and then click on the left side on *SESSIONS RESULTS AVERAGE FOR DIFFICULTY LEVEL 30*.

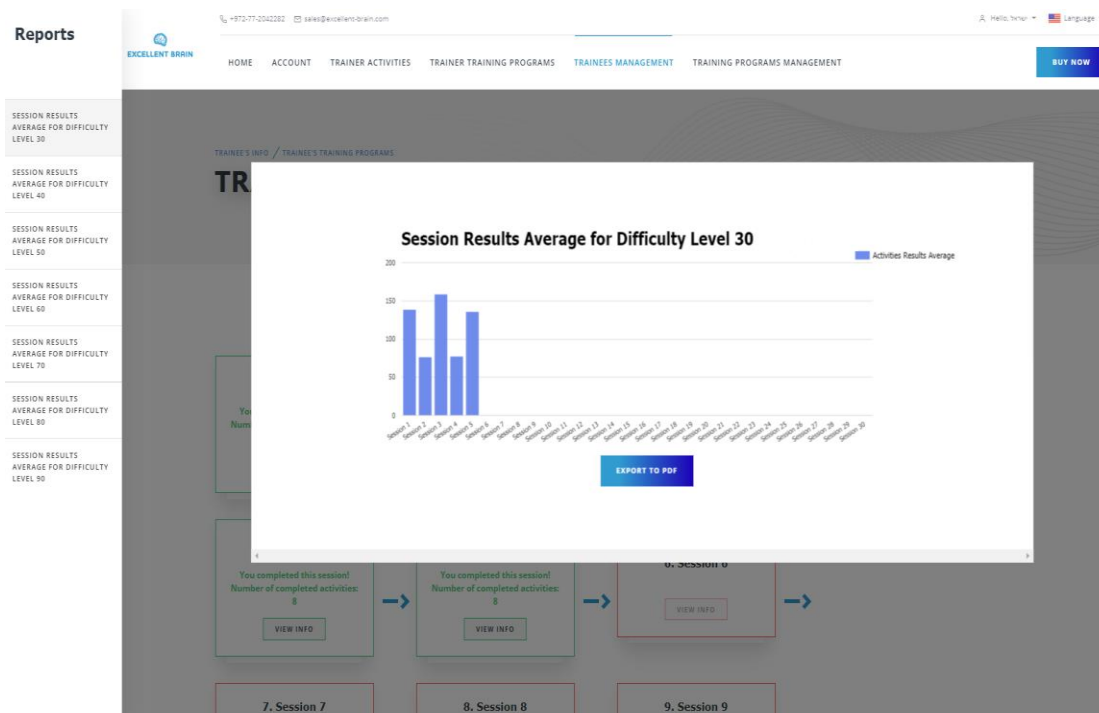

11. **View Session Results** – In order to view a specific Session results and notes, reach the page at section 9 and then choose the required Session by clicking on *VIEW INFO*.

Reports

EXCELLENT BRAIN

HOME ACCOUNT TRAINER ACTIVITIES TRAINER TRAINING PROGRAMS TRAINEES MANAGEMENT TRAINING PROGRAMS MANAGEMENT BUY NOW

ACTIVITIES SCORES

TRAINER'S TRAINING PROGRAMS / TRAINEE'S SESSIONS

SESSION 1 - MICHAL MICHAELI

#	Activity Name	Difficulty Level	Protocol	Score in this session	Best score for this level of activity	Date	Notes	
1	Cupid	30%	Attention	117	117	11/07/2019	<input type="text"/>	SAVE
2	My Earth	30%	Attention	160	160	05/06/2019	<input type="text"/>	SAVE
3	Shark	40%	Attention	142	142	05/06/2019	<input type="text"/>	SAVE
4	Launch the Rocket	40%	Attention	288	288	05/06/2019	<input type="text"/>	SAVE
5	Resize the Diamond	50%	Attention	153	153	05/06/2019	<input type="text"/>	SAVE
6	Lift the Brain	50%	Attention	104	104	11/07/2019	<input type="text"/>	SAVE
7	Fireworks	50%	Attention	68	68	05/06/2019	<input type="text"/>	SAVE
8	Glowing Cube	40%	Attention	64	64	05/06/2019	<input type="text"/>	SAVE


Total score: 1096

SAVE SESSION NOTES

12. View a Report Displaying Session Activities Results Compared to the Same Activities in the Same Difficulty Level from Previous Sessions - Reach the page at section 11 and then click on the left side bar on *ACTIVITIES RESULTS*.

You will get this graph showing in red the current result of this activity in the difficulty level, and in blue the maximum result of the same activity in the same difficulty level from previous Sessions.

Since this is only the first Session, they are equal.

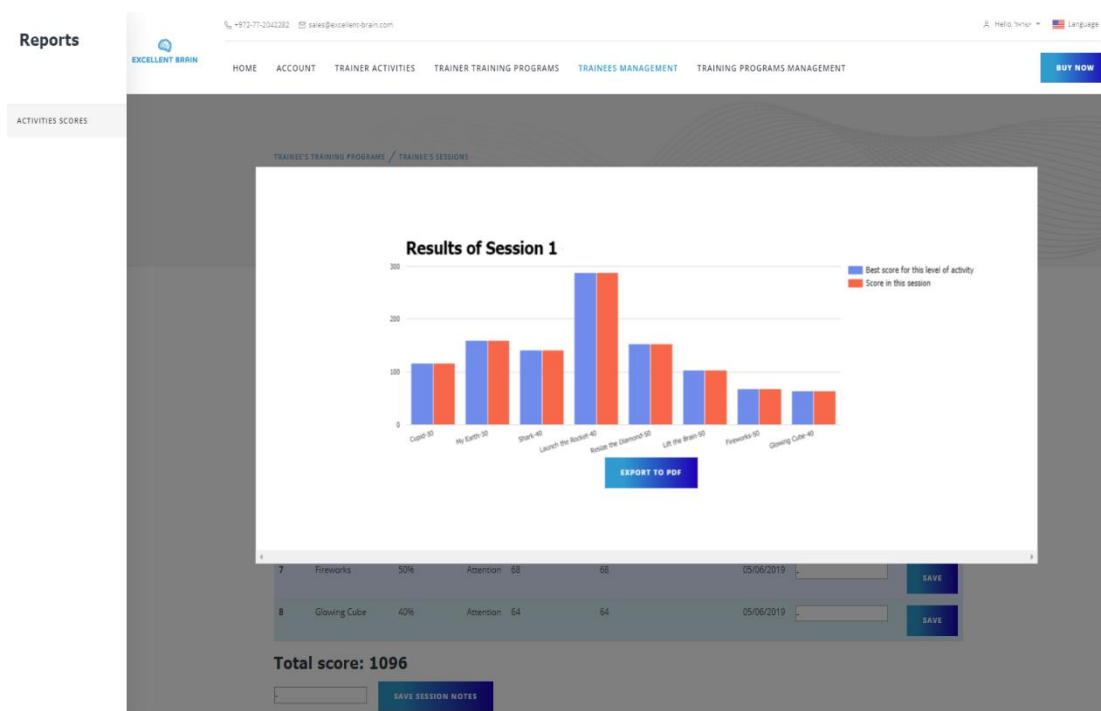

13. **View Your Personal Training Programs Management** – In order to create personalized Training Programs for trainees, click at the top of the page on *TRAINING PROGRAMS MANAGEMENT* and the you will be moved to a page with this content –

- Your Personal Sessions List (Top Section) – You can create up to 10 Sessions for each trainee and another 10 Sessions for yourself.
- Your Personal Training Programs List (Middle Section) – You can create up to 1 Training Program for each trainee and another one for yourself.
- Your Trainees List (Bottom Section)

a. Begin from Creating New Sessions.

At the **top section** named *YOUR SESSIONS* click on the green button named *CREATE NEW SESSION*. A window will open where you need to fill in the following fields-

- Session Name
- Session Description
- Number of Session Activities
- Activities Names and their Difficulty Levels (30-90)

Click Finish and then you can see that the Session was successfully added to the table.

Important Information:

- In order to view the details of the Session you created, click on the blue arrow named *Show Activities* under Actions on the right side of the required Session. A list of Activities of this Session will open.
- In order to edit the details of the Session you created, click on the yellow pencil named *Edit* under Actions on the right side of the required Session. A window similar to the one we created the Session from will open.
- In order to completely delete the Session you created, click on the red bin named *Delete* under Actions on the right side of the required Session. A window will open where you will be required to approve the deletion. This action is irreversible.

You can add as many Sessions until reaching the quota (see above in this page, 10 Sessions for each trainee and another 10 Sessions for yourself).

b. Continue to Creating New Training Programs.

At the **middle section** named *YOUR TRAINING PROGRAMS* click on the green button named *CREATE NEW TRAINING PROGRAM*. A window will open where you need to fill in the following fields-

- Training Program Name
- Training Program Description
- Number of Sessions
- List of Sessions Names to Add

Click Finish and then you can see that the Training Program was successfully added to the table.

Important Information:

- In order to view the details of the Training Program you created, click on the blue arrow named *Show Sessions* under Actions on the right side of the required Training Program. A list of Sessions of this Training Program will open.
- In order to edit the details of the Training Program you created, click on the yellow pencil named *Edit* under Actions on the right side of the required

Training Program. A window similar to the one we created the Training Program from will open.

- In order to completely delete the Training Program you created, click on the red bin named *Delete* under Actions on the right side of the required Training Program. A window will open where you will be required to approve the deletion. This action is irreversible.

You can add as many Training Programs until reaching the quota (see above in this page, 1 Training Program for each trainee and another one for yourself).

At this point all of these Training Programs will be available to you on Excellent Brain Software (version 0.95.10 and above).

c. Continue to Assigning Training Programs to Trainees.

At the **bottom section** named *YOUR TRAINEES*, mark the required trainee at the left side of the table and click on the green button named *ASSIGN TRAINING PROGRAMS TO ALL CHOSEN TRAINEES*. A window will open where you need to choose which Training Programs will be assigned to the chosen trainee. Click Finish and then the process is complete. This action can also be done by clicking the yellow plus sign named *TRAINING PROGRAM ASSIGNING* under Actions, on the right side of the chosen trainee. A window will open where you need to choose which Training Programs will be assigned to the chosen trainee.

Important Information:

- In order to view the assigned Training Programs of a certain trainee, click on the blue arrow named *Show Training Programs* under Actions on the right side of the required trainee. A list of the trainee's assigned Training Programs will open.

- In order to completely remove all assigned Training Programs of a certain trainee, click on the red bin named *Remove All Training Programs* under Actions on the right side of the required trainee. A window will open where you will be required to approve the removal. This action is irreversible.

TRAINING PROGRAMS MANAGEMENT

CREATE NEW SESSION		DELETE		Your Sessions		
				Allowed Sessions quota 2 / 80		
<input type="checkbox"/>	#	Session Name	Description	Number of Activities	Add/Update Date	Actions
<input type="checkbox"/>	1	Trainer Session 1	Description for Trainer Session 1	2	11/07/2019	+ - X
X						
#	Activity Name		Difficulty (%)			
1	Fireworks		30%			
2	The Golden Mandela		40%			
<input type="checkbox"/>	2	Trainer Session 2	Description for Trainer Session 2	3	11/07/2019	+ - X
Showing 2 out of 2 Sessions						
				Previous	1	Next

CREATE NEW TRAINING SESSION		DELETE		Your Training Programs		
				Allowed Training Programs quota 1 / 8		
<input type="checkbox"/>	#	Training Program Name	Description	Number of Sessions	Add/Update Date	Actions
<input type="checkbox"/>	1	Trainer Training Program 1	Description for Trainer Training Program 1	2	11/07/2019	+ - X
X						
#	Session Name	Description	Number of Activities	Add/Update Date		
1	Trainer Session 1	Description for Trainer Session 1	2	11/07/2019		
2	Trainer Session 2	Description for Trainer Session 2	3	11/07/2019		